

# ASSI® Hidalgo Lower Blepharoplasty Retractor



The Hidalgo Blepharoplasty Retractor facilitates lower eyelid fat transportation. It provides superior visibility for release of the *arcus marginalis* and exposure of the orbital rim. It enables accurate needle passage deep to the *orbicularis* muscle when placing percutaneous fixation sutures.




Non-Insulated Retractor  
ASSI.AG97326  
6.5cm Regular Model

ASSI.AG98126  
14.5cm Long Model

Insulated Retractor  
ASSI.AG97626  
6.5cm Regular Model

ASSI.AG98226  
14.5cm Long Model

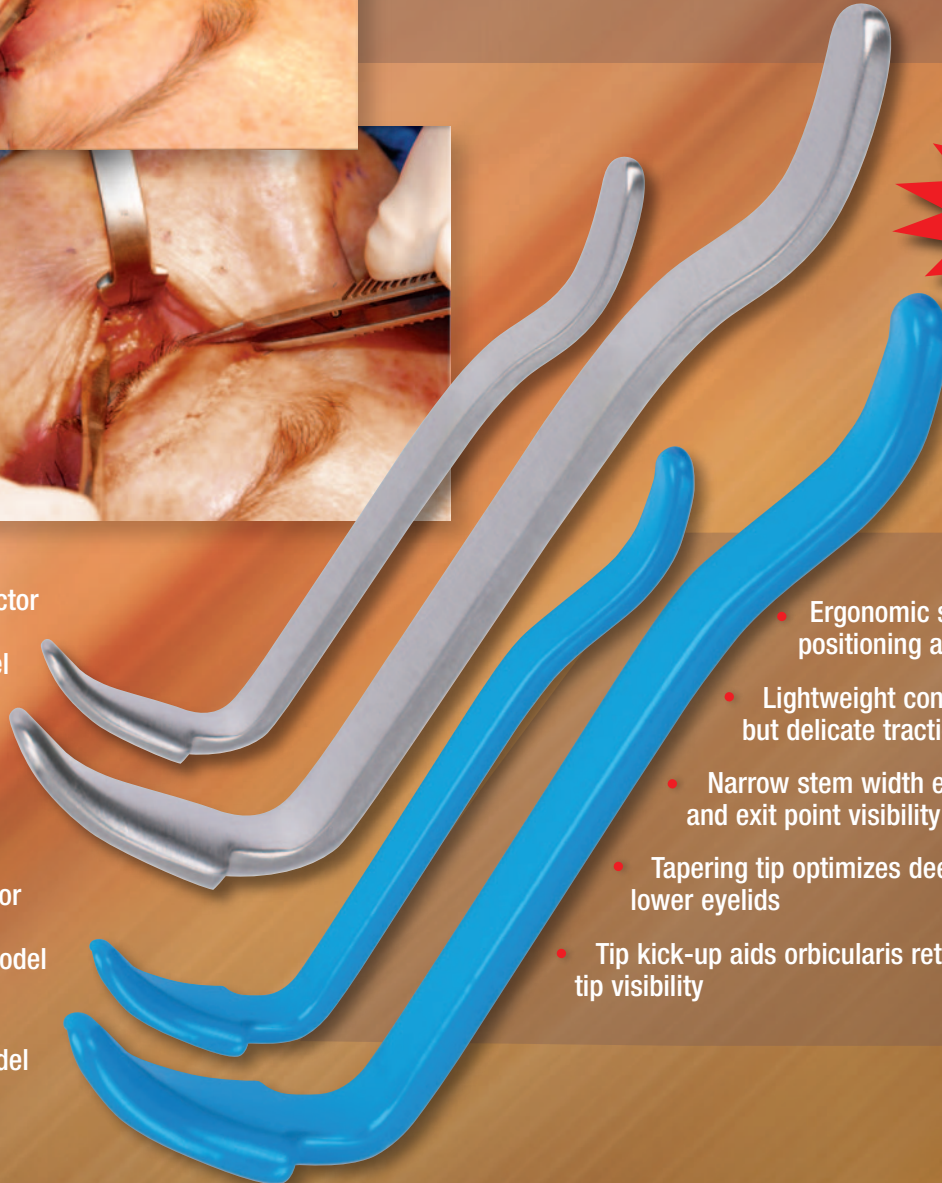

- Ergonomic stem design facilitates positioning and stability

- Lightweight construction permits firm but delicate traction force application

- Narrow stem width enhances suture entrance and exit point visibility

- Tapering tip optimizes deep exposure in tight lower eyelids

- Tip kick-up aids orbicularis retraction and needle tip visibility

Designed by:  
David A. Hidalgo, M.D.  
Clinical Professor of Surgery  
Cornell Medical Center, NY, NY

**assi**®   
ACCURATE SURGICAL & SCIENTIFIC INSTRUMENTS®

accurate surgical & scientific instruments corporation

300 Shames Drive, Westbury, NY 11590

800.645.3569 516.333.2570 fax: 516.997.4948 west coast: 800.255.9378

Info: [assi@accuratesurgical.com](mailto:assi@accuratesurgical.com) • Orders: [orders@accuratesurgical.com](mailto:orders@accuratesurgical.com) • [www.accuratesurgical.com](http://www.accuratesurgical.com)

Not all ASSI products shown in our literature or on our website are available for sale in Canada